



The **Wheel**® is a uniquely effective device. It is physician designed to give immediate temporary relief from pain and soreness of muscle related disorders.

## Physiology

Strength, speed, flexibility and endurance are known qualities of healthy muscle. A muscle's ability to prolong these desirable qualities is directly dependent on an unimpeded, nutrient-rich and oxygen-laden blood supply.

The **Wheel**® provides the many documented benefits of trigger point therapy and myofascial release. Its purpose is to extinguish myofascial trigger points and augment the qualities of healthy muscle.

## Indications

The **Wheel**® is successfully used in the musculotendinous therapy of:


- carpal tunnel syndrome • occipital neuralgia
- tension headache • TMJ disorder
- epicondylitis • tendinitis
- and a host of other myofascial lesions.

CPT code 97140 • E 0935



## Experience the Magic



- Hold the handlegrip in the palm of your hand...apply pressure with the thumb and index finger.
- The **Wheel**® is rolled lengthwise over the involved muscle, using short back & forth passes.
- Use the **Wheel**® to search and find myofascial trigger points...these areas are tight, knotty, ropy and tender.
- After locating the trigger point, continue the rolling process for about 30 seconds.
- On a pressure-pain scale of 10 – use 5...the procedure should produce a "feel-so-good" pain.
- Chronic, stubborn trigger points may initially require several short, therapeutic sessions per day.
- Your goal is a healthy muscle free of pain, tightness, tenderness and other signs of offending trigger points.

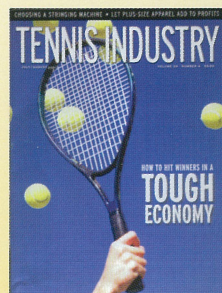

**"The Wheel can effectively warm up muscles before workouts."**



**Relief is just a Roll-A-Way!**

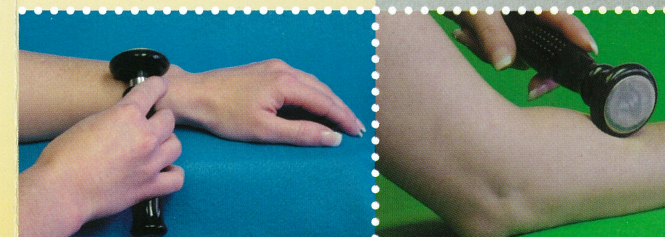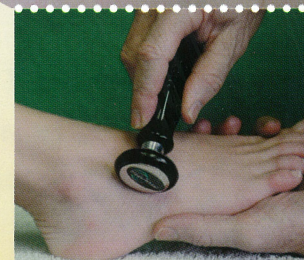





***"It works!  
Ask my husband, Robbie."***

Genie Windholz • Atlanta, GA

***"The Wheel allows me to spend  
more hours writing and editing  
than I otherwise could."***

John C Lowe

The Fibromyalgia Research Foundation  
Houston, TX

***"I recommend the Wheel for  
all my patients with symptoms of  
carpal tunnel syndrome."***

Sydney Wade, PT • Athens, Greece

***"I dispense the Wheel  
for my patients to use  
at the job site."***

R G Hotchkiss, DC • Rockford, IL



Manufactured by

**RPI** of Atlanta

Atlanta, GA USA  
[www.intracell.net](http://www.intracell.net)  
Patents Granted and Pending

ICT 04/304

